

DIGITAL X-RAY SYSTEM

ERAY SMART FLOOR MOUNTED 10~80kW



WIDE APPLICATION IN

Clinics | Polyclinics | Nursing
Homes | ICU | OT | Sports
Centers | Trauma Centers



ENHANCE YOUR WORKFLOW & EFFICIENCY

General-Purpose Radiography System

The ERAY RAD Series is a state-of-the-art digital X-ray system that offers high-quality diagnostic imaging with minimal radiation exposure. Its compact and space-saving design allows for easy installation in small spaces, making it a perfect choice for clinics and private practices with limited space. Additionally, its patient-friendly design provides a comfortable environment for patients during examinations, contributing to an overall positive experience. With its advanced features and user-friendly interface, the ERAY RAD Series is a reliable and efficient digital X-ray system suitable for various clinical settings.

SOFT HANDLING


ERAY SMART offers soft handling. Allowing the operator to move the tube head unit (THU) with minimal physical pressure, reducing fatigue.

CONSISTENTLY CLEAR IMAGES

ERAY SMART provides dependable and excellent imaging for patients with diverse anatomical structures, ensuring consistently clear results with more than 1280 APR.

VERTICAL CHEST BUCKY STAND

The floor-mounted chest bucky stand is a versatile and flexible system, designed for easy maneuvering and handling. Equipped with a 17"x17" dynamic flat panel detector, it can be used for lay down, upright, or oblique x-ray examinations. The chest bucky stand features a manually operated tilting design, which allows for easy adjustment in 6 steps of 15 degrees, catering to a wide range of radiography examination.

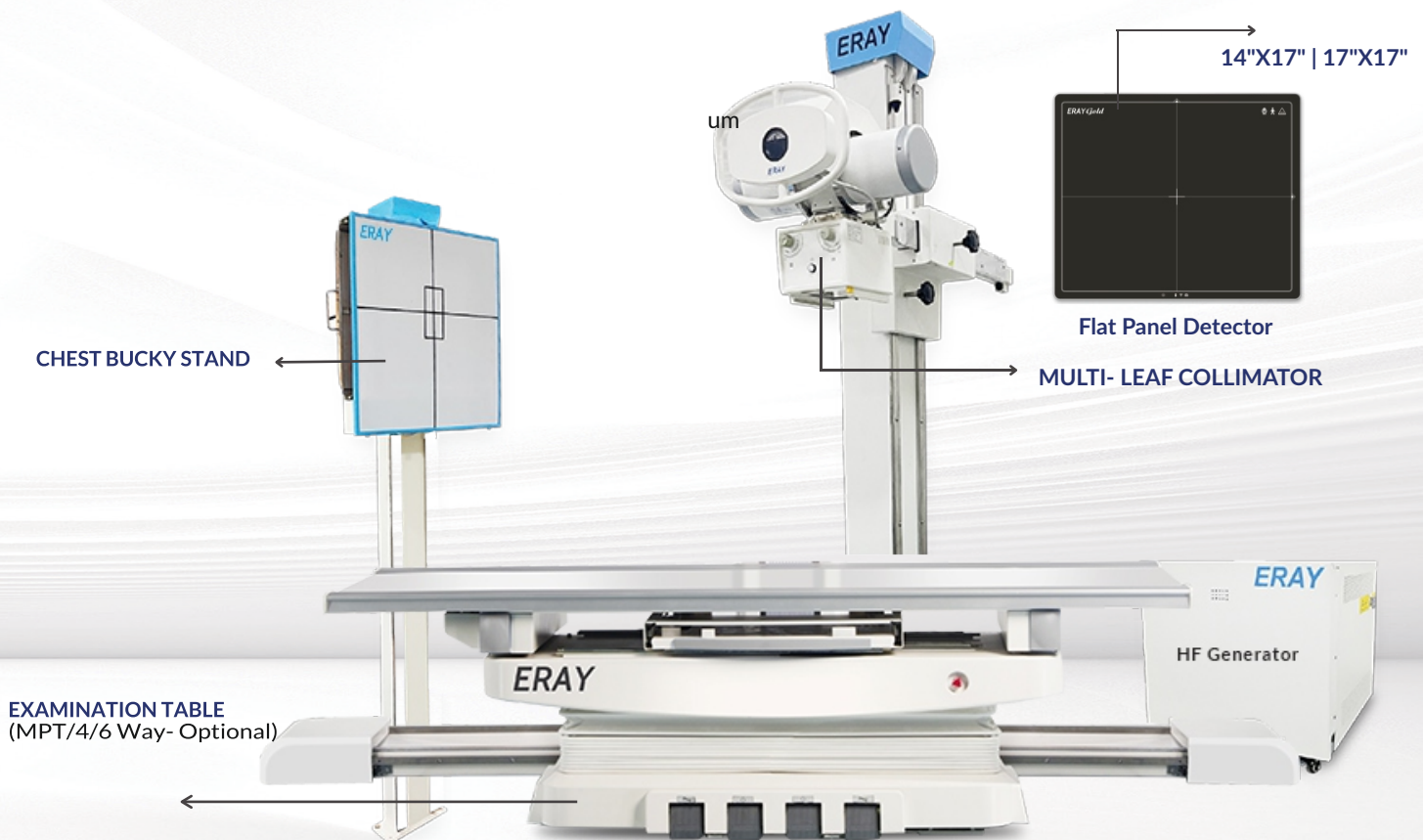


STANDARD FEATURES

- Digital display exposure parameters.
- Parameter increase / decrease switches
- Ready and X-Ray ON switch with indicators.
- Bucky selection switch
- A dual action hand switch with retractable cord for radiation protection of operator.

ERAY GOLD FPD

- Resolution 3072X3072
- 140um
- 70cm Drop Height
- Weight 3.2
- Waterproof IP56

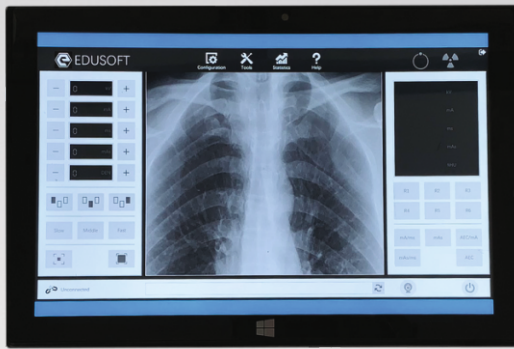

X-RAY CONSOLE

The system includes a 10" touch screen tablet that enables the operator to precisely position the 4-point display (kV, mA, mAs & Time) with ease. This feature allows the operator to focus on patient care, as the tablet can be easily reached.




THE HF ADVANTAGE

- ✓ Lower Radiation Doses
- ✓ Diagnostic Precision
- ✓ High Quality Images
- ✓ Less Skin Dose

GLOBAL NETWORK

Edusoft Proudly Exporting to Many Countries Like GCC Bangladesh, Dubai, Nepal, Singapore, Indonesia, Malaysia, + Latin America, Philippines, Africa

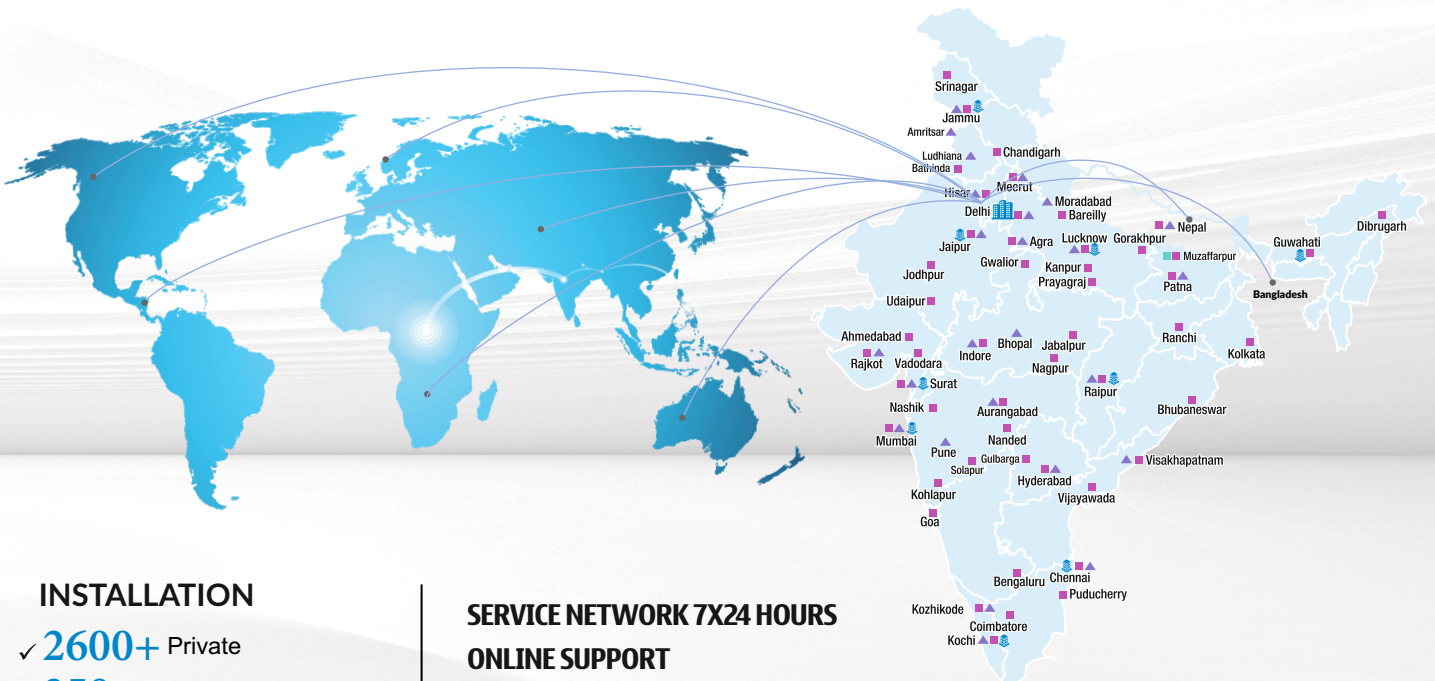

INSTALLATION

- ✓ **2600+** Private
- ✓ **350+** Government

50+ Channel Partners / Dealers

37+ Field Service Engineers

50+ Sales + Mkts

SERVICE NETWORK 7X24 HOURS ONLINE SUPPORT



Professional Training



Remote Support



Regular Visits



HEAD OFFICE

DEALERS

SERVICE ENGINEERS

Technical Highlights

ERAY X-RAY

MODEL	ERAY RAD 200S	ERAY RAD 250S	ERAY RAD 300S	ERAY RAD 500S	ERAY RAD 630S	ERAY RAD 1000S
GENERATOR (kW)	10kW	20kW	32kW	40kW	50kW	80kW
INVERTER FREQUENCY	220khz	220khz	220khz	220khz	220khz	220khz
LINE NOMINAL, PHASE	Three Phase, 380V, 50/60Hz	Three Phase, 380V, 50/60Hz	Three Phase, 380V, 50/60Hz	Three Phase, 380V, 50/60Hz	Three Phase, 380V, 50/60Hz	Three Phase, 380V, 50/60Hz
LINE VOLTAGE RANGE	±10% 50/60 Hz					
kV RANGE	40kV- 125kV	40kV- 125kV	40kV- 125kV	40kV- 125kV	40kV~150kV	40kV~150kV
mA RANGE	10mA-200mA	10mA-250mA	10mA-300mA	10mA-500mA	10mA-630mA	10mA-1000mA
TIMER RANGE	0.01sec-10sec					
mAs RANGE	0.1-200mAs	0.1-250mAs	0.1-320mAs	0.1-320mAs	0.1-320mAs	0.1-320mAs
MAX. POWER OUTPUT	125kVp~80mA 100kVp~100mA	125kVp~160mA 100kVp~200mA	125kVp~250mA 100kVp~320mA	150kVp~260mA 80kVp~500mA	150kVp~333mA 80kVp~630mA	150kVp~660mA 80kVp~1000mA



X-RAY TUBE

FOCAL SPOT SIZE	1.0/2.0 mm	1.0/2.0 mm	1.0/2.0 mm	1.0/2.0 mm	0.6/1.5mm	0.6/1.2mm
MAX. ANODE kHU	140kHU	140kHU	140kHU	140kHU	200kHU	300kHU 400kHU(Optional)
MAX.kV	125kV	125kV	125kV	125kV	150kV	150kV

ACCESSORIES

X-RAY TABLE (Standard)	Horizontal Bucky Table
X-RAY TABLE (Optional)	Multi-Position
WALL VERTICAL BUCKY	Tilting (0° to 90° & -20°)
CHEST BUCKY STAND	Optional
X-RAY ACCESSORIES	Optional

Clinical Images

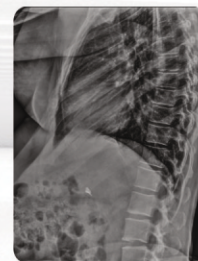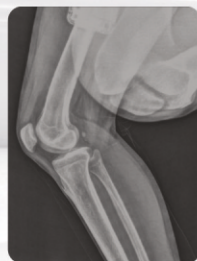

NOTE: *For reference purpose only

Edusoft Healthcare Pvt. Ltd.

📍 B-83, Phase-II, Mangolpuri Industrial Area, New Delhi-110034, India

✉️ ceo@edusofthealth.com

☎️ 1800-120-280-280

📞 +91 8851248073, +91 8826066621

• Noida • Lucknow • Chandigarh • Pune • Bhopal • Chennai • Bengaluru • Kolkata
• Hyderabad • Jaipur • Guwahati • Ahmedabad • Mumbai

Copyright © 2025 Edusoft Healthcare Pvt. Ltd. All rights reserved. Specifications and designs are subject to change without notice. Non-metric weights and measurements are approximate. All data were deemed correct at time of creation. Edusoft is not liable for errors or omissions. All brand, product, service names and logos are trademarks and/or registered trademarks of their respective owners and are hereby recognized and acknowledged.



Scan Code or visit
www.edusofthealth.com
to learn more